



Custom measurement form: PCO ORTHOSES CARD

PATIENT ID


PRODUCT CODE

NOTES



REH4MAT
MANUFACTURER OF MODERN ORTHOPAEDIC BRACES
36-060 Głogów Małopolski,
ul. Piaski 47, POLAND
Adres korespondencyjny
ul. Truskawkowa 17, Widna Góra
37-500 Jarosław, POLAND
tel. +48 (016) 621 42 20
biuro@reh4mat.com

MEDICAL
DEVICE
CLASS 1
WYRÓB
MEDYCZNY
KLASA 1

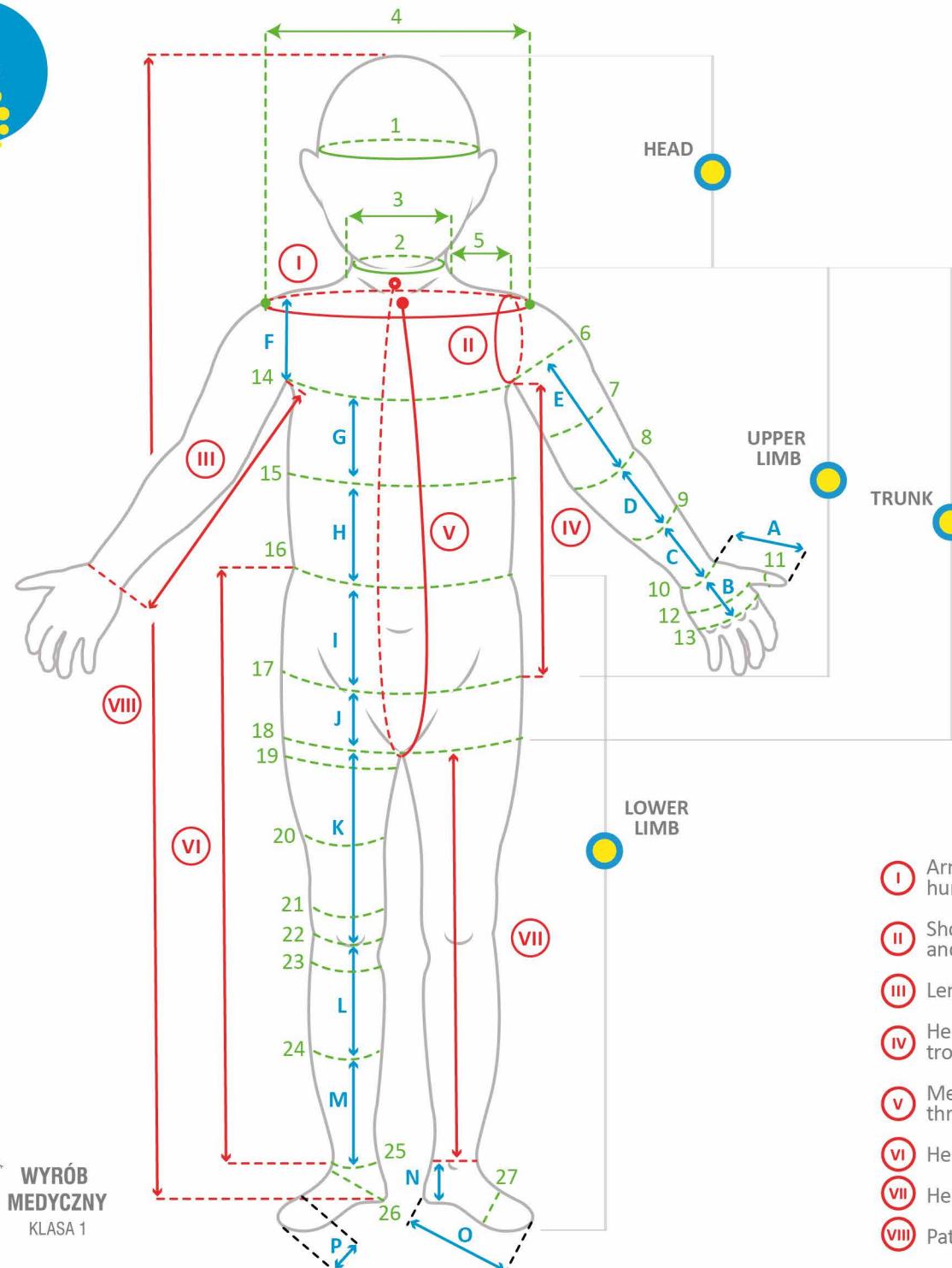


WARNING!

- Chest and upper limb measurements should be taken with patient lying down.
- Hip and lower limb measurements are taken with patient standing.
- All measurements should be taken in centimeters!
- If the finger put between the child's body and the PCO product meets resistance, the size of the product has been selected correctly.

IMPORTANT!

Only a correctly filled in measurement card is a guarantee of proper compression of the product. We are not responsible for wrongly taken measurements!

HOW TO MEASURE?

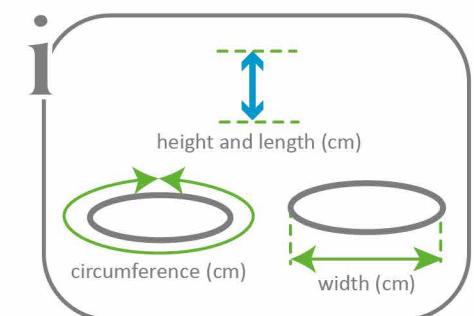

- I Arms circumference at the level of the head of the humerus
- II Shoulder circumference measured through the armpit and acromion
- III Length measured from armpit to wrist
- IV Height measured from the armpit to the greater trochanter of the femur
- V Measurement from the manubrium of the sternum through the perineum to the cervical vertebra C7
- VI Height measured from waist to ankle
- VII Height measured from perineum to ankle
- VIII Patient is height

REQUIRED MEASUREMENTS

PRODUCT	MEASUREMENT	LENGTH MEASUREMENT	CONTROL MEASUREMENT	ZIP
HEAD				
PCO-H-01	1			
TRUNK				
PCO-T-01	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, V, VIII	
PCO-T-02	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, V, VIII	
PCO-T-03	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, V, VIII	
PCO-T-04	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, V, VIII	
PCO-T-09	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, VIII	V
PCO-T-10	2, 3, 4, 5, 6, 7, 8, 14, 15, 16, 17, 18	E, F, G, H, I, J	I, II, III, IV, VIII	V
PCO-T-11	2, 3, 4, 5, 6, 7, 8, 9, 10, 14, 15, 16, 17, 18	C, D, E, F, G, H, I, J	I, II, III, IV, VIII	V
PCO-T-16	2, 3, 4, 5, 6, 14, 15, 16, 17, 18	F, G, H, I, J	I, II, IV, VIII	V
PCO-T-17	2, 3, 4, 5, 6, 7, 8, 14, 15, 16, 17, 18	E, F, G, H, I, J	I, II, III, IV, VIII	V
PCO-T-18	2, 3, 4, 5, 6, 7, 8, 9, 10, 14, 15, 16, 17, 18	C, D, E, F, G, H, I, J	I, II, III, IV, VIII	V
PCO-T-23	2, 3, 4, 5, 6, 14, 15, 16, 17, 18, 19	F, G, H, I, J	I, II, IV, V, VIII	V
PCO-T-24	2, 3, 4, 5, 6, 7, 8, 14, 15, 16, 17, 18, 19	E, F, G, H, I, J	I, II, III, IV, V, VIII	V
PCO-T-25	2, 3, 4, 5, 6, 7, 8, 9, 10, 14, 15, 16, 17, 18, 19	C, D, E, F, G, H, I, J	I, II, III, IV, V, VIII	V
UPPER LIMB				
PCO-A-01	6, 7, 8, 9, 10, 11, 12, 13	A, B, C, D, E	III	
PCO-A-02	8, 9, 10, 11, 12, 13	A, B, C, D	III	
PCO-A-03	9, 10, 11, 12, 13	A, B, C	III	
PCO-A-07	6, 7, 8, 9, 10	C, D, E	III	
PCO-A-08	8, 9, 10	C, D	III	
LOWER LIMB				
PCO-L-01	16, 17, 18, 19, 20, 21, 22, 23, 24, 25	I, J, K, L, M	VI, VII, VIII	
PCO-L-02	16, 17, 18, 19, 20, 21, 22, 23, 24	I, J, K, L	VI, VII, VIII	
PCO-L-03	16, 17, 18, 19, 20, 21	I, J, K	VI, VII, VIII	
PCO-L-08	19, 20, 21, 22, 23, 24, 25	K, L, M	VII	
PCO-L-09	19, 20, 21, 22, 23, 24	K, L	VII	
PCO-L-10	19, 20, 21	K	VII	
PCO-L-11	23, 24, 25	L, M	VII	
PCO-L-12	24, 25, 26, 27	M, N, O, P	VII	

No.	MEASUREMENT	CIRCUMFERENCE	WIDTH	LENGTH
1	HEAD CIRCUMFERENCE	(cm)	(cm)	(cm)
2	NECK CIRCUMFERENCE			
3	NECK WIDTH			
4	SHOULDERS WIDTH			
5	SHOULDER-NECK WIDTH			
6	ARM CIRCUMFERENCE			
7	CIRCUMFERENCE BETWEEN 6 AND 8			
8	ELBOW CIRCUMFERENCE WITH A SLIGHTLY BENT HAND			
9	CIRCUMFERENCE BETWEEN 8 AND 10			
10	WRIST CIRCUMFERENCE			
11	THUMB CIRCUMFERENCE			
12	HAND CIRCUMFERENCE			
13	MCP JOINTS CIRCUMFERENCE			
14	MEASUREMENT UNDER THE ARMPITS			
15	MEASUREMENT BETWEEN 12 AND 14			
16	WAIST CIRCUMFERENCE			
17	HIP CIRCUMFERENCE			
18	MEASUREMENT AT THE PERINEUM HEIGHT			
19	GROIN CIRCUMFERENCE			
20	CIRCUMFERENCE BETWEEN 19 AND 20			
21	CIRCUMFERENCE 5 CM ABOVE THE KNEECAP			
22	KNEE CIRCUMFERENCE			
23	CIRCUMFERENCE 5 CM UNDER THE KNEECAP			
24	CALF CIRCUMFERENCE BETWEEN 22 AND 25			
25	ANKLE CIRCUMFERENCE			
26	HEEL CIRCUMFERENCE			
27	METATARSUS CIRCUMFERENCE			

LEGEND:  - MEASUREMENT UNAVAILABLE

CONTROL MEASUREMENT	
I	
II	
III	
IV	
V	
VI	
VII	
VIII	
ZIP	
	FRONT
	BACK